

HIFU VAGINAL TIGHTENING MACHINE




3.0mm / 4.5mm

FEMALE PRIVATE CARE

- » PRIVATE TIGHTEN
- » PRIVATE LUBRICATION
- » PRIVATE SENSITIVE
- » SENSITIVE HEALTH


THEORY

The Private Tightening HIFU System uses a noninvasive ultrasonic focusing technique to directly focus on the mucosal lamina and muscle layer. Using ultrasonic waves as the energy source and taking advantage of its penetration and focusability, the system will send out ultrasonic energy focusing in the lamina propria and muscle fiber layer. In a predetermined depth layer, a higher intensity of ultrasonic region, called focus region, is formed. In 0.1 second, the temperature of the region can reach to above 65 °C, so the collagen is reorganized in a predetermined depth layer and the normal tissue outside the focal region is undamaged. To ensure the desired depth layer tissue can obtain the ideal effect of collagen contraction, reorganization and regeneration. Ultimately, the effect of vagina tightening is achieved.

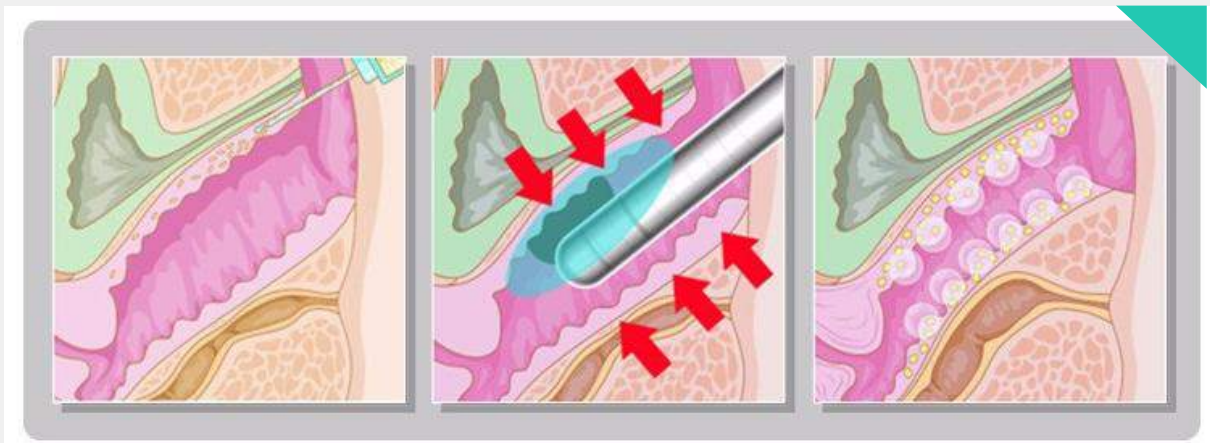

APPLICATION

1. Potent firming effect to improve sexual life
2. Private lubrication of increasing secretion and eliminating dryness
3. Improving private sensitivity by cell regeneration for reaction enhancement
4. Comprehensive improvement of private health and reducing infection

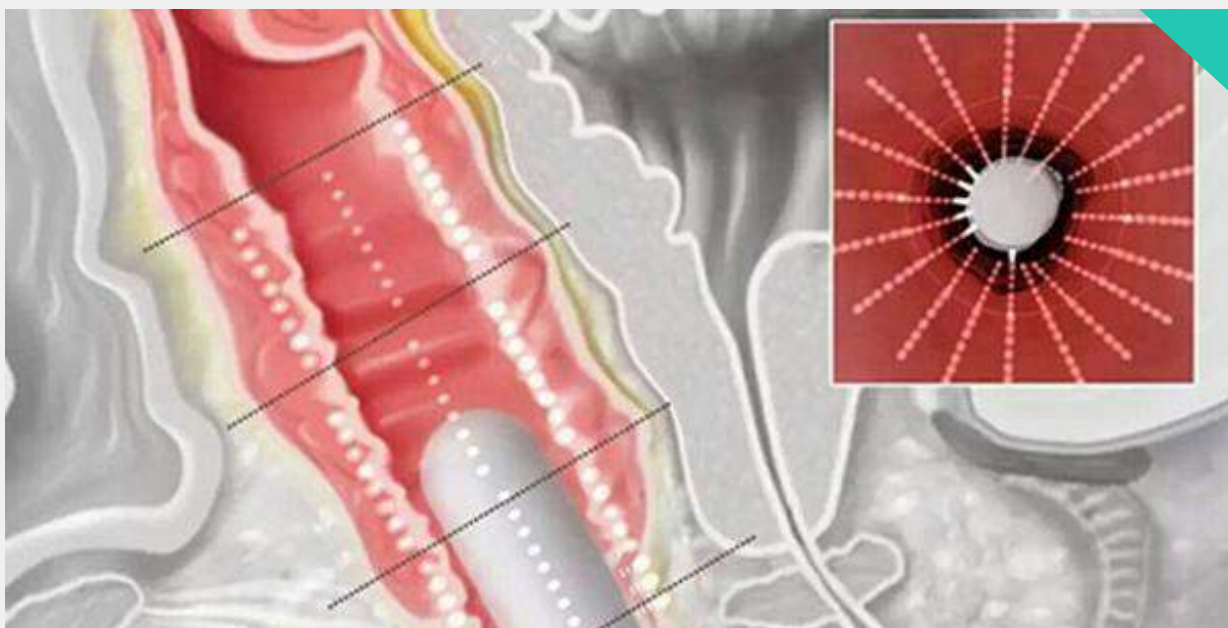

HANDLE AND EMITTER





Handle

3.0mm Emitter

4.5mm Emitter

Soft Baffle

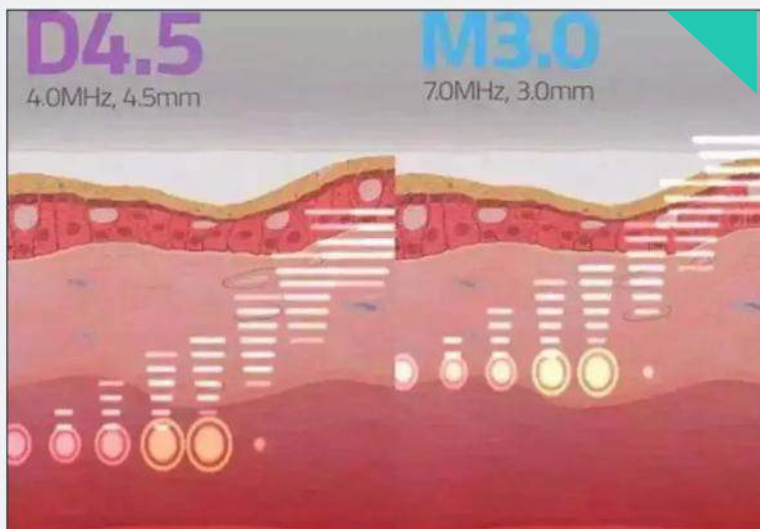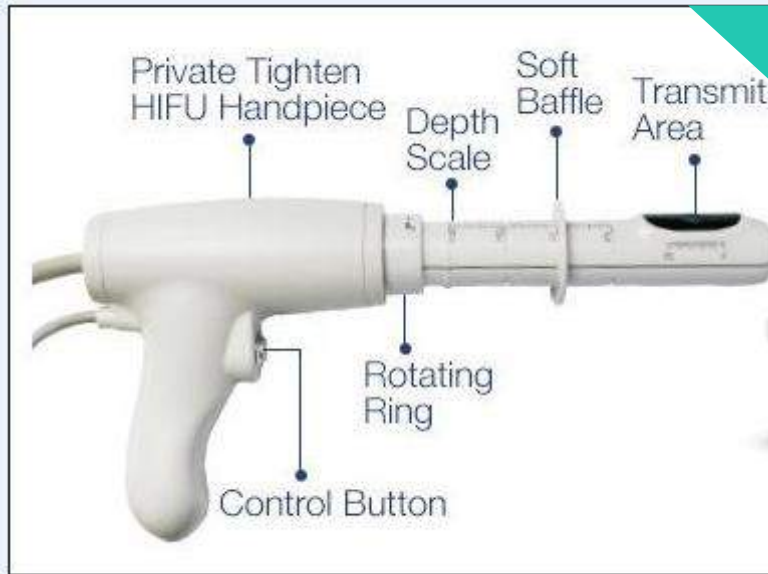

DS 3.0mm (7MHz):

Muscle Layer

DS 4.5mm (4MHz):

SMAS

ADVANTAGES

- 1.The equipment has two different working heads for different operational effect: 3.0mm and 4.5mm.
- 2.It has precise depth adjustment system.
- 3.360°rotation emission: total care for vaginal.
- 4.Precise depth adjustment system.
- 5.Non-invasive, no downtime, no recovery time, can have sexy in 3 days after treatment.
- 6.It is easy and convenient for operation, and no consumable items are required, which greatly saves the treatment cost.
- 7.Tightening and shaping effect is obvious after treatment. It can be maintained at least 18 to 24 months after one treatment and realize negative growth of skin age once a year.
- 8.The normal life and work will not be influenced if you make up the moment.
- 9.Easy and convenient: 20 minutes treatment can get instantly firmer, less treatment times, simple operation.

SPECIFICATIONS

Energy Type	HIFU (high intensity focused ultrasound)
Emitter	DS-3.0mm, DS-4.5mm
Power	0.2-2.0J
Spacing	1.0-10mm
Length	5-25mm
Rotating Angle	5-25°
Voltage	AC100V-240V, 50/60Hz

CONTACT INFORMATION

Get in touch

www.letemps.paris

Tel: +33 783 5883 68



3 years free warranty, lifetime maintenance,
24 hours online service

Machine will come with
.User's manual
.Maintenance manual
.Training DVD